55" 4K Medical Display

FM-A5504DS FM-A5504DG









- Excellent color and detail reproduction on a large screen.
- Wide selection of signal input choices.
- Built in multi-window layouts.



See Every Detail

This monitor delivers ultra-sharp picture clarity with maximum viewing angles, allowing

several clinicians to view high levels of detail without loss of image quality or distortion. Its 3,840 x 2,160 resolution greatly enhances images from endoscopic camera systems. Input source options include DVI, HDMI, SDI, and DisplayPort. The FM-A5504DG model adds 12G SDI to the impressive selection of signal choices.






Multipl

Multiple Windows

There is plenty of room on the FM-A5504D monitor series for multi-window requirements. Choose from PIP, PBP, Triple, or Quad layouts that are preset and ready-to-go in the on-screen display menus.

Edge-to-edge protective front glass gives this monitor a slim, sleek appearance, and allows for easy cleaning of the front surface. Rounded bezel corners help reduce bumping injuries, especially in today's active operating room environments. A back cable cover is provided with each monitor to manage cords and wires.

FM-A5504DS ON/OFF SDI 3 DVI 4 SDI 1 DVI 2 USB DISPLAY DISPLAY PORT 1 5V AC INLET SWITCH IN IN IN IN OUT 2 RS232C (SVC) PORT 2 IN OUT • •) • • **GROUND** SDI 4 DVI 3 SDI 2 DVI 1 GPIO HDMI 5V OUT 1 IN IN IN

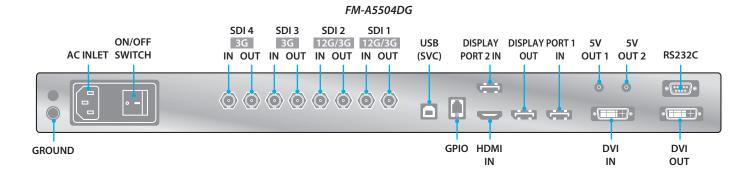

55" 4K Medical Display

FM-A5504DS FM-A5504DG




Specifications


General Features

Item		Description
Panel		55 inch TFT LCD (LED)
Resolution		3840 x 2160 pixels
Aspect Ratio		16:9
Active Area		1209.60 (H)mm x 680.40 (V)mm
Pixel Pitch (mm)		0.315 x 0.315
Response Time (typical)		8 ms (gray to gray)
Number of Colors		1.07 Billion
Brightness (typical)		500 cd/m ²
Gamut		BT.709 and BT.2020 compatible
Contrast Ratio (typical)		1100:1
Surface Treatment		Anti-glare
Viewing Angle (CR>10)		R/L 178°, U/D 178°
Input Signal	FM-A5504DS	1 x HDMI 2.0 2 x DP 1.2 (SST) 4 (Quad) x DVI (single link) 4 (Quad) x SDI (3G)
	FM-A5504DG	1 x HDMI 2.0 2 x DP 1.2 (SST) 1 x DVI (single link) 4 (Quad) x SDI (3G), 2 x SDI (12G)


FSN Medical Technologies products meet EU RoHS Directive 2011/65/ EU+2015/863/EU and Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH) obligations. FSN manages the information flow as required by regulations.

Item		Description
Output Signal	FM-A5504DS	1 x DP 1.2 (SST)
	FM-A5504DG	1 x DP 1.2 (SST) 1 x DVI (single link) 4 (Quad) x SDI (3G), 2 x SDI (12G)
Power Supply		SMPS (AC 100 ~240V)
Power Consumption		140W max
Accessories Included		User manual, AC power cord, DP cable, HDMI cable, DVI-D cable, BNC cable, Remote control and batteries
Unit Dimension		1292.5(W) x 770.5(H) x 93.5(D) mm 50.89(W) x 30.33(H) x 3.68(D) inch
Package Dimension		1450(W) x 930(H) x 305(D) mm 57.09(W) x 36.61(H) x 12(D) inch
IP Rating		IP30 - overall, IP42 - front side
IK Rating		IK07 with 5H AR glass
\\\ : .	FM-A5504DS	36.65 kg, 80.80 lbs. (monitor with cover) 49.73 kg, 109.64 lbs. (shipping package)
Weight	FM-A5504DG	36.65 kg, 80.80 lbs. (monitor with cover) 50.23 kg, 110.74 lbs. (shipping package)
Certifications		ANSI/AAMI ES60601-1 (UL60601-1), CAN/CSA C22.2 No. 60601-1, FCC Class B, EN60601-1, EN60601-1-2, CE, MDR(EU) 2017/745 Class I Medical Device

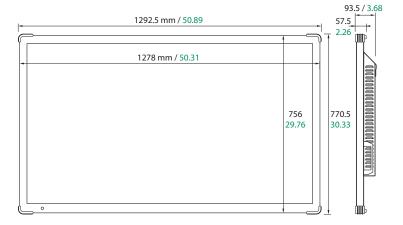

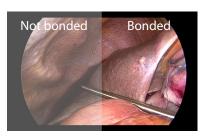

FSN's optical bonding process (right) improves contrast ratio and brightness, and increases humidity resistance and breakage resistance.

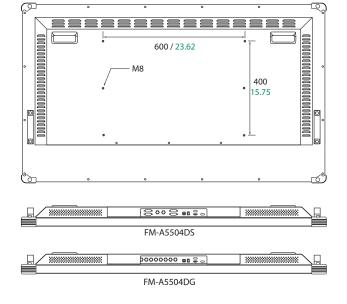

Specifications are subject to change with or without notice. Doc. # FSN2032 Rev. 11/21

www.fsnmed.com

Foreseeson Custom Displays, Inc. 2210 E. Winston Road Anaheim, CA 92806 USA Tel: 714-300-0540 Fax: 714-300-0546

Aventura, FL 33180 USA Tel: 714-507-3855

Foreseeson Korea 404B, Pangyoinnovalley B 253 Pangyo-ro Bundang-gu Seongnam-si Gyeonggi-do, Korea 463-400 Tel: +82-31-8018-0780

Fax: +82-31-8018-0786

Foreseeson GmbH Industriestrasse 38a 63150 Heusenstamm, Germany Tel: +49 6104 64398 0 Fax: +49 6104 64398 11

Foreseeson UK Ltd. 1 Wolsey Road East Molesy Surrey KT8 9EL **United Kingdom** Tel: +44 (0) 208 546 1047 Foreseeson (Shanghai) Medical Equipment Co., Ltd. Room 307, 3F No. 56, 461 Hongcao Road Caohejing Development District Xuhui, Shanghai 200233 Tel: 86-21-6113-4188